

Protecting Patient Privacy, Data Integrity and Caregivers' Time



Secure access to Connex® vital signs devices with Imprivata®

Clinicians need fast access to patient vitals. But hospitals must secure medical devices to protect **patient privacy, data integrity** and **compliance**.

That's why the Welch Allyn Connex® family of vital signs devices feature secure access with Imprivata.

The process is simple.

Clinicians open a session with Imprivata, either at the EMR workstation or elsewhere. This authentication is extended to the medical device, and access is granted in a single step to take patient vitals.




Scan or tap badge at bedside—no username or password required. (RFID or barcode scanning options available).



Record vitals as an authorised caregiver.

Caregiver ID Data Path

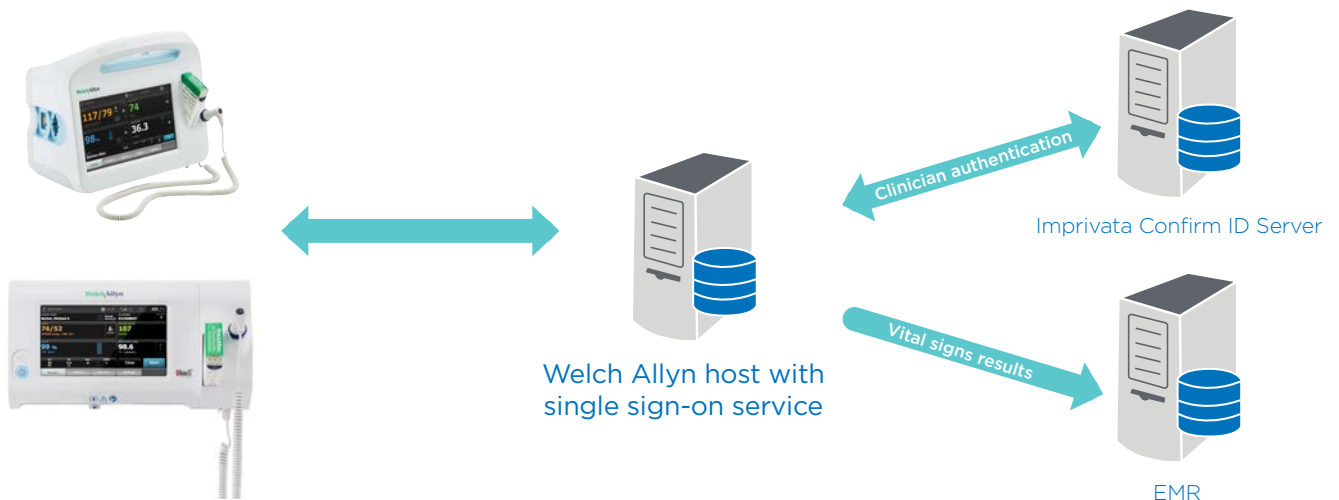

Protect What Matters.

Extending your badge-access method from desktops to Connex devices can help you:



Improve data integrity into the EMR.

Many hospitals are increasing protections to ensure only authorised users can add information to the EMR. The Welch Allyn - Imprivata integration enables fast, secure authentication for accessing and transacting patient information—ensuring the user is authorised to update the patient chart.



Control access to the device.

A simple badge scan or tap enables quick authentication at the bedside, granting device access only to authorised users—and helping protect patient health information.



Streamline clinical workflows.

Keep the security of a username and password without the frequent entry. After signing in through Imprivata in accordance with hospital policies, caregivers simply scan or tap their badge to log into Connex devices—so they can spend their valuable time on patient care.

It's all about protecting what matters to you—patient privacy, data integrity and caregivers' time.

Caregiver ID Type Supported	Connex Spot Monitor	Connex Vital Signs Monitor
Barcode scanning	●	●
RFID (13.56 MHz, unencrypted cards)	●*	

Welch Allyn Host Solutions Supported
Cerner VitalsLink
Welch Allyn Network Connectivity Engine 5.1X and higher
Welch Allyn Connex CS (CSAS) 1.8.1 and higher

*RFID support on Connex Spot Monitor software version 1.32 and higher

Contact your Welch Allyn representative or visit www.welchallyn.co.uk to learn more.

Welch AllynUK Ltd.
 Clinitron House, Excelsior Road,
 Ashby-de-la-Zouch
 LE65 1JG
 0207 365 6780
ukcustomer care@welchallyn.co.uk

WWW.WELCHALLYN.CO.UK



WelchAllyn®

Hill-Rom®

©2018 WELCH ALLYN, INC. ALL RIGHTS RESERVED. DOC. NO: SM4184UK REV.B MCI5478, AUG 2018

The Connex Vital Signs Monitor and the Connex Spot Monitor are intended to be used by clinicians and medically qualified personnel for monitoring vital signs of all patients. Medical device (93/42/EEC); Class IIa Notified Body: DQS Medizinprodukte GmbH (CE0297) Manufacturer: Welch Allyn, Inc., 4341 State Street Road, Skaneateles Falls, NY 13153 USA. EC Rep.: Welch Allyn Limited, Navan Business Park, Dublin Road, Navan, County Meath, Republic of Ireland. These medical devices are a regulated health product which, pursuant to such regulation bears a CE mark. Welch Allyn recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices. Welch Allyn reserves the right to make changes without notice in design, specifications and models. The only warranty Welch Allyn makes is the express written warranty extended on the sale of its products.